

Go for comfort.



swissmade 

The most advanced automatic reusable pen.

- Spring driven delivery for an effortless and comfortable injection process
- Intuitive and proven twist motion dose dialling and push-button activation for maximum patient acceptance
- Easy dose setting and correction ensuring high patient comfort
- Bayonet connection and sliding lead screw for quick cartridge change
- For insulin and other multidose therapies
- High quality look & feel including metal outer housing




More confidence. More success. With Ypsomed Delivery Systems.




ServoPen® – automatic reusable pen

Intuitive spring driven delivery


The ServoPen® is an intuitive spring driven insulin pen providing an effortless, yet familiar handling experience.

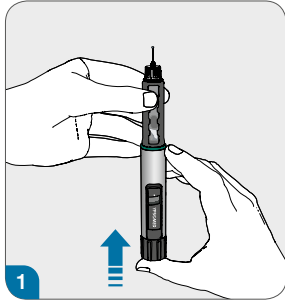

1 Prime

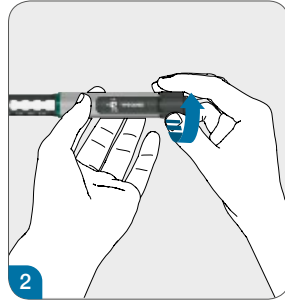

2 Dial dose

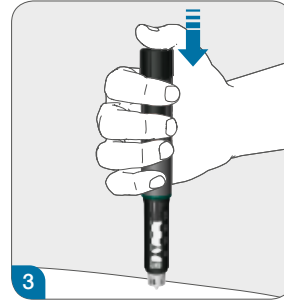

3 Inject



Advanced pen technology

Innovative automatic delivery based on a proven handling concept

Handling

- The shortest available stroke distance on the market for a comfortable delivery
- Enlarged dose display for good readability for visually impaired patients
- Cartridge holder with bayonet connection and sliding lead screw makes changing the cartridge quick and easy
- Audible and tactile feedback during dosing, dose correction and delivery to maximise patient confidence
- Real-time display providing the possibility to interrupt the injection and to know how much still has to be delivered
- Pre-injection remaining dose indication ensuring correct dosing up to the last unit

Customization and manufacturing

- Off the shelf device easily customized for your drug
- The ServoPen® can be customised for different fill volumes and dose ranges
- Robust design to minimise complaint rates and allow a pen lifetime of 3 years



Dose dialling by twisting is intuitive and easy. Enlarged dose display provides easy to read numbers.



Design and customisation

Your drug deserves a unique presentation

The ServoPen® can be customised with different colours for the cap and outer housing and with customer specific printing.




The ServoPen® can be supplied with exchangeable colour rings to colour-code pens when a patient uses different drugs in different ServoPens®.



Dimensions and specifications

Facts and figures

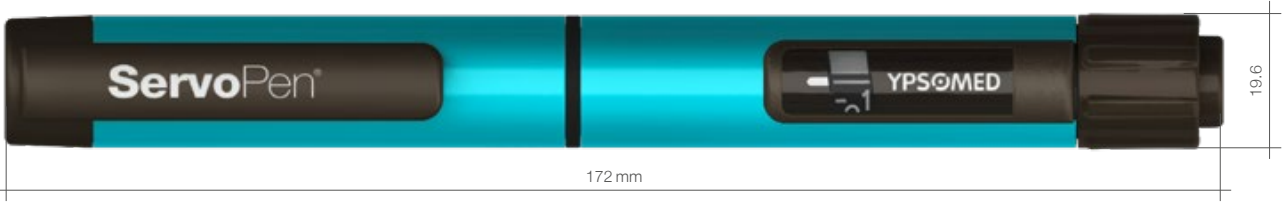


Device specifications

Primary container	Standard 3 ml cartridges according to ISO 11608-3*	Min. dose/ dose increment	0.01 ml
Fill volume	3 ml	Max. dose	0.60 ml
Number of dose increments	60	Compatible needle type	Compatible with standard and safety pen needles** – mylife™ Clickfine® and mylife™ Clickfine® AutoProtect™ pen needles recommended

* 3 ml cartridge inner diameter of 9.7 mm recommended
 ** Pen needles according to ISO 11608-2

Cartridge dimensions (mm)

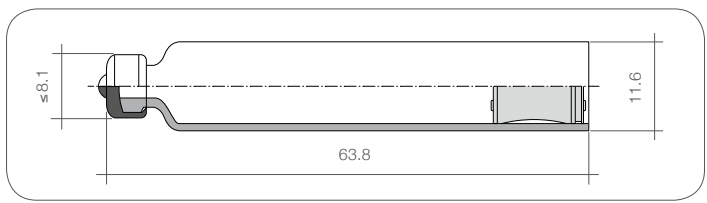

Your delivery system

Drawing from our experience and competence

- Ypsomed has extensive experience in injection device development and manufacturing of more than 30 years and has numerous devices on the market.
- You can expect high quality support for choosing the primary container, customising of the industrial design, demonstrating dose accuracy and preparing the relevant device documentation for the registration of your product.

More **confidence.**
More **success.**
With **Ypsomed Delivery Systems.**



Pen systems



Autoinjector systems



Patch injector systems



Smart services

Ypsomed Delivery Systems provides a complete range of drug delivery products and services to biopharmaceutical companies. We offer everything from development and design to manufacturing and packaging, giving patients and customers more confidence and leading to more market success.